



homecare equipment

Kirton Postural Care Seating

Florien II



homecare equipment

For the semi and non-ambulant user including side transfers
Mid to high posture and pressure care needs

A dynamic and highly mobile seating system.

The Florien II enables **greater independence** by promoting **correct posture** and facilitating a **supportive** comfortable seated position.

The **removable arm** is a fantastic feature of this chair making it suitable for those who are hoisted or able to side transfer.

Tilt-in-Space positioning - mechlock operated

Legrest elevation - mechlock operated

Removable arms - easy side transfers

Moulded backrest - enhanced comfort and seating experience


Integral head support - articulated

25° negative legrest angle:

- for standing transfers or to accommodate knee contractures

Three seat length options - 435mm | 460mm | 485mm

Weight limit -180kg

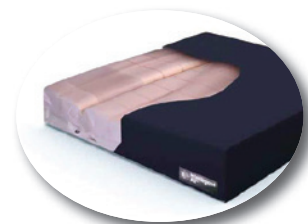








Includes Intelligent Air Cushion



Adjustable Lateral Supports



Optional Flip-up Footrest

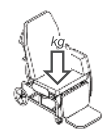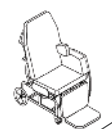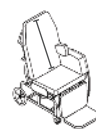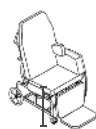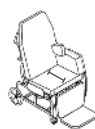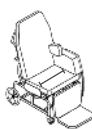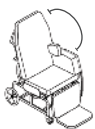

Tilt-in-space	Seat Widths	Seat Depths	Seat Height	Legrest Elevation	Back Height	Overall Width	Weight Capacity
5° - 35°	400mm 450mm 500mm	435mm 460mm 485mm	550mm	15° - 90°	830mm	700mm 740mm (Large Wheel)	180kg

Florien Elite



homecare equipment

For the ambulant, semi and non-ambulant user
Mid to high posture and pressure care needs

With both tilt-in-space and back angle recline functionality, the Florien Elite has a choice of a **waterfall** or **contoured** backrest and **three armrest options** to meet the needs of different care settings.

Dynamic upholstered arms float across the seat to achieve the appropriate width for each user. The arm can also swivel backwards for side transfers or easier hoist access.



Dynamic Upholstered
Basic Cushioned



Dynamic Drop Armrest
Swivel Function



- Manual tilt-in-space - Independent back recline
- Manual elevating legrest - negative retraction
- Multi-positioning headrest - integral
- Adjustable seat - width | depth | height
- Contoured or Waterfall backrest
- Armrest style options
- Basic cushioned armrest for
 - portering | drip-stand | phlebotomy
- Reflexion Foam™ - Dartex on contact surfaces
- Weight limit - 160kg



Tilt-in-space	Seat Widths	Seat Depths	Seat Height	Legrest Elevation	Back Height	Overall Width	Weight Capacity
40°	400 - 500mm (12mm increments)	430 - 530mm	550mm 575mm 600mm	-10° - 70°	800mm	668mm	160kg

G - 2




homecare equipment

For the semi and non-ambulant user
High posture and pressure care needs

Underpinned by **scientific research** and ergonomic design principles, the G-2 is a groundbreaking model which achieves new levels of **comfort** and **immersion**.

This elegant chair has been designed to ensure all key adjustments are **simple** and **effective**, meeting the needs of both users and carers. Perfect for those seated for lengthy periods and who require support and positioning assistance.

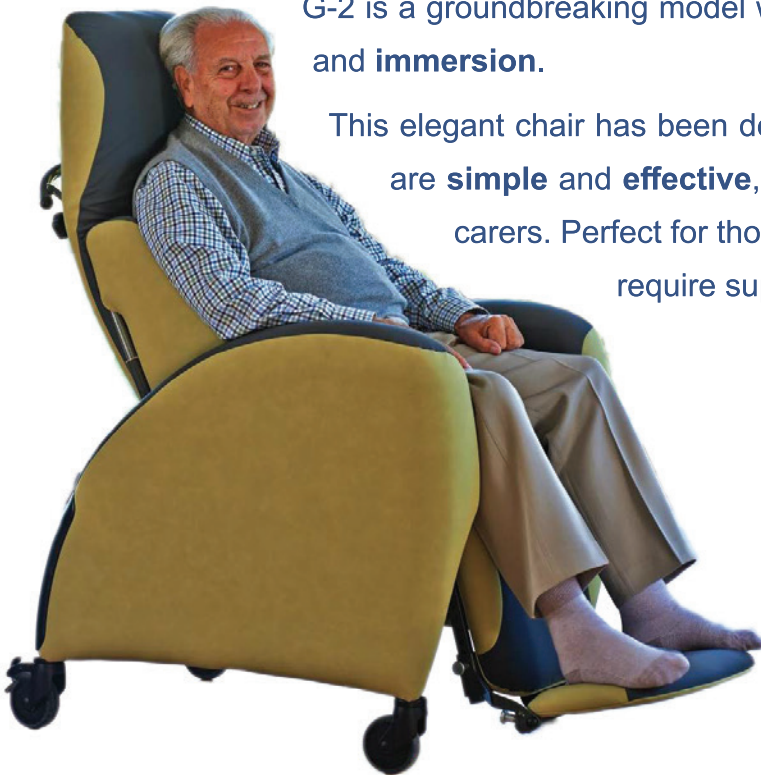

Regulated Motion

(optional with powered version)

Automated tilt repositioning between 10° - 25° ensuring dynamic movement and continuous changing of pressure hotspots



Tilt-in-Space and Legrest Elevation - manual | electrically powered

Articulating supportive backrest

Helical motion lateral supports

Arcing legrest motion - maintains footboard positioning


Adjustable seat length: 380mm - 480mm

Adjustable footboard height and angle

Includes Dyna-Tek Intelligent air cushion

Optional leg-laterals

Weight limit - 180kg

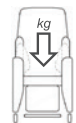


Tilt-in-space	Seat Widths	Seat Depths	Seat Height	Footboard Height	Back Height	Overall Width	Weight Capacity
50°	400 450mm 500 550mm	380 - 480mm	530mm	380 - 540mm	800mm	700mm 800mm	180kg

For the semi and non-ambulant user including side transfers
Mid to high posture and pressure care needs

The **Advanced Comfort Backrest**, incorporating Reflexion™ Foam, moulds itself to the contours of the user, helping **disperse body weight** and providing **superior comfort** and support. A **Waterfall Back** is also available for more cushioned, immersive support.


The chair can be **raised or lowered** to attain the correct **seat height** for each user and helps maintain an **even weight distribution** by ensuring feet are firmly positioned on the floor. This also assists the ambulant user to stand from the chair.

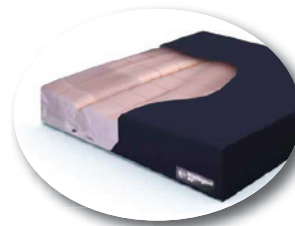

Includes
Intelligent Air
Cushion



- Tilt-in-Space positioning - gas action | electrically powered
- Back angle recline - gas action | electrically powered
- Legrest elevation - gas action | electrically powered
- Includes Dyna-Tek Intelligent air cushion
- Adjustable seat - height and length
- Advanced Comfort Backrest - Reflexion™ Foam on contact surfaces
- Integral headrest - multi-positioning
- Integral lateral supports - adjustable
- Weight limit - 160kg

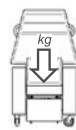

Tilt-in-space	Seat Widths (Manufactured to spec)	Seat Depths	Seat Height	Legrest Elevation	Back Height	Overall Width	Weight Capacity
26° - Manual 30° - Electric	400 450mm 500 550mm	420 - 540mm	430 - 550mm	0° - 90°	850mm	700mm 800mm	160kg



For the semi and non-ambulant user
For those with Huntington’s disease and other **neurological conditions**



Originally designed for those with **Huntington’s disease** and other **neurological conditions**, the Omega is **comfortable** and **secure**, yet sufficiently **robust** to withstand severe movements while providing a **safe** and **relaxed** seated position.


The sculpted back and contoured seat provides good positioning for activities such as eating and drinking, while relaxing and sleeping is encouraged by the comfort of the uniquely linked back and seat recline - as the back reclines, the seat angle increases.

- Manual back recline and variable seat angle
- Pivoting seat front edge - retain knee to floor height
- Robust, durable frame - modular back | arm | seat components
- Horizontally and vertically sculpted back - aids midline stability
- Two back height options: 740mm | 800mm
- Three seat width and length options
- Includes Intelligent Air cushion
- Supportive Reflexion™ Foam legrest - removable
- Weight limit - 150kg



Optional
Sliding Footrest

Reflexion™ Foam
Removable Legrest

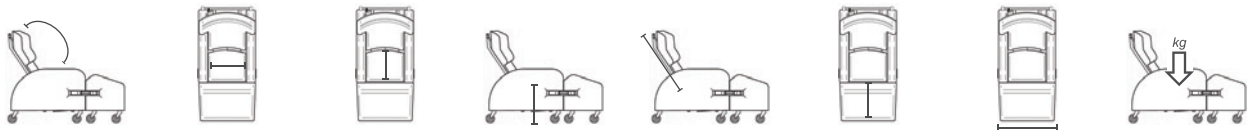

Recline Tilt	Seat Widths	Seat Depths	Seat Height	Back Height	Legrest Depth	Overall Width	Weight Capacity
32° - 57° (Back) 20° - 35° (Seat)	380mm 430mm 510mm	430mm 470mm 500mm	485mm	740mm 800mm	460mm	670mm 750mm	150kg

For the non-ambulant user

Profound posture and pressure care needs



Dynamic features provide the ultimate in **positioning** and **postural management** for individuals who are at risk of developing an **unconventional body shape**, unable to maintain a **midline position** and require **additional support**. A versatile seating system with adjustability to ensure the core principles of seating are met for each individual. Available in **three chair sizes** to ensure optimum fit.



Regulated Motion

(optional with powered version, not available with Duo Mini)

Automated tilt repositioning between 10° - 25° ensuring dynamic movement and continuous changing of pressure hotspots



Tilt-in-Space and Legrest Elevation - manual | electrically powered

Adjustable seat - length | width | height* (*not on manual version)

3 Chair sizes available - Major | Minor | Mini (paediatric)

Simple tool-free adjustments - back height | footboard | arm height


Arcing legrest motion - maintains flip-up footboard positioning

Removable arms for side transfer access

Split-wing shoulder support - adjustable head & shoulder support

Adjustable headrest - precise head positioning

Two back style options - waterfall | sculpted

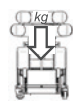

	Tilt-in-space	Back Recline	Negative Legrest	Seat Widths	Seat Depths	Seat Height	Back Height	Overall Width	Weight Capacity
Major	35° (28° - manual version)	90° - 130°	-25°	350 400mm 450 500mm	380 - 500mm	*510 - 630mm (*manual version)	670 - 940mm	680 - 830mm	160kg
Minor	35° (28° - manual version)	90° - 130°	-25°	300 350mm 400 450mm	320 - 460mm	*500 - 620mm (*manual version)	610 - 820mm	620 - 770mm	160kg
Mini	40°	95° - 130°	n/a	200 250mm 300 350mm	240 - 340mm	430 - 530mm	490 - 660mm	570 - 640mm	57kg



You can use the following **SELECTION MATRIX** of features to help you find the right dynamic chair... though a full assessment is always recommended.

Some features are common to all the chairs in the range:

Modular design - refurbish and replace over time

Postural support and pressure relief

Simple to use | Easy to clean

Versatility - domestic use | multi-user environments | clinical settings | care settings

	Florien II	Florien Elite	G - 2	Encora	Omega	Duo
Removable side	✓	✓				
Footrest height adjustment	✓	✓	✓	✓		✓
Legrest elevation	✓	✓	✓	✓		✓
Integral pressure relief	✓	✓	✓	✓	✓	✓
Tilt-in-space - Manual	✓	✓	✓	✓		✓
Tilt-in-space - Electric			✓	✓		✓
Regulated Motion			✓			✓
Back rest recline		✓		✓		✓
Seat height adjustment		✓		✓		Electric Only
Seat depth adjustment	Depth options available	✓	✓	✓	Depth options available	✓
Seat width adjustment	Width options available	✓	Width options available	Width options available	Width options available	✓
Back height adjustment					At time of manufacture	✓
Lateral Support	Optional	Optional	Built-in	Built-in	Built-in	Optional
Waterfall back option	✓	✓		✓		✓
Weight capacity	180kg	160kg	180kg	160kg	150kg	160kg - Major 160kg - Minor 57kg - Mini

* Other fabrics available but this will result in longer delivery times and possible increased costs

*To ensure the pressure relieving properties of your chair are working to their full potential, we prefer Dartex on all user **contact surfaces**.*

Dartex is a strong, multi stretch, vapour permeable, waterproof fabric which is simple to wipe clean and truly complements the pressure relief incorporated into the chairs.

