



homecare equipment





ProCair Quad

Delivering the highest level of comfort with 1 in 4 air-cell alternation as well as superior safety features including a fully sealed base, the ProCair Quad is especially suitable for hospital and aged care clients with a high to very high risk of pressure injuries.



ProCair Quad

AcuteCair

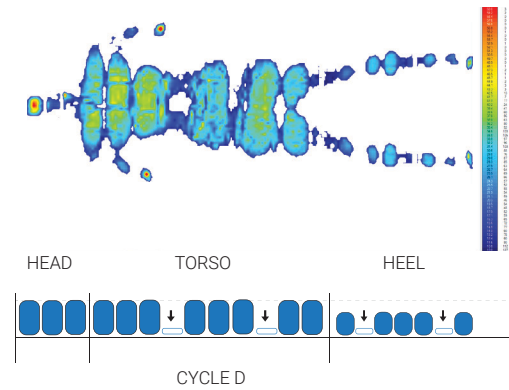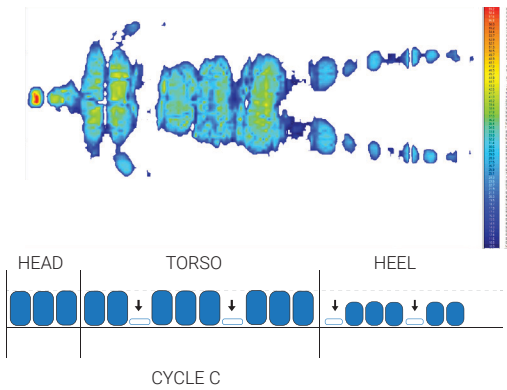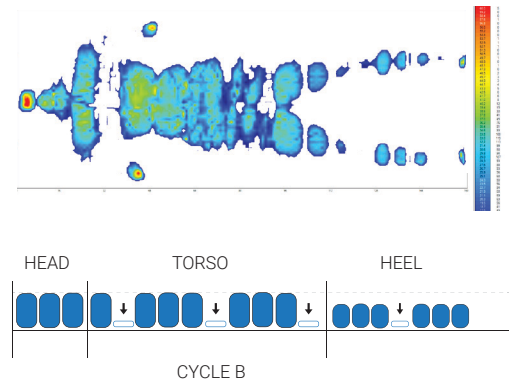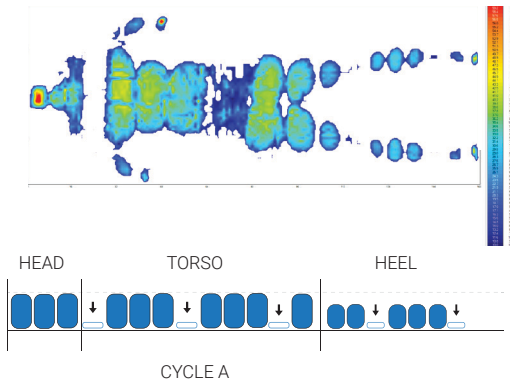

1-in-4 Alternation


- Supports 75% patient's body during active therapy
- Redistributes the pressure and delivers greater comfort

Cell in cell technology

- Lower air cells remain inflated during therapy which maintains support of patient, provides superior infection control, and prevents patients from bottoming out

Comfort Control

- 3 Comfort levels to suit individual need

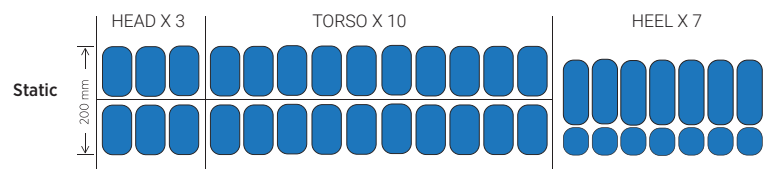

ProCair Quad

Features

Multi-stretch water resistant, vapour permeable cover

To protect against shear and friction and reduce moisture against skin

Bed exit alarm

Audible alarm can be activated to manage falls risk

TPU Cells

Automatically adjust to the ambient temperature of the room

Heel Zone

- Individually detachable heel zone cells to allow complete off-loading of the vulnerable calcaneus
- Tapered heel sloped design with narrower cells

AutoCair

Automatic weight sensing ensures correct weight setting for effective therapy

Backup battery

Allows for continuous active therapy up to 3 hours in the event of a power failure and during patient transport

24 Hours transport mode

Rapid CPR

Release in case of emergency

Fully sealed base

For maximum infection control



Automatic fowler boost

Protects patient from bottoming out when positioned with the head raised between 15-45 degrees



Cable management system

Ensure cords are neatly secured along either side of the mattress to reduce risk of trips and potential damage to power cords



Side Bolsters

Provide stability for the patient during sit-to-stand transfers and bed mobility



homecare equipment

ProCair Quad

Specifications

	SINGLE	KING SINGLE
CODE	APMPC-R04	APMPC-R04K
CAPACITY (SAFE WORKING LOAD)	40-250kg	40-250kg
SYSTEM ALTERNATION	1 in 4	1 in 4
WARRANTY	3 Year	3 Year
PRESSURE INJURY RISK	High to Very High	High to Very High
COMPLIANCE	IEC60601-1, IEC60601-1-2 and IEC60601-1-11	IEC60601-1, IEC60601-1-2 and IEC60601-1-11
ARTG	289458	289458
MATTRESS DIMENSIONS (L x W x H)	2000 x 880 x 200mm	2000 x 1050 x 200mm
MATTRESS WEIGHT	12.5kg	14kg

Refer to the Technical Data Sheet for detailed specifications.

Accessory



Fast Inflation/deflation pump

Setup or deflate mattress in under 2 minutes

